

SL350 - SL450

ALL THE ESSENTIALS FOR EYE EXAM



SL350 - SL450



SLIT LAMPS

NO AKL: KEMENKES RI AKL 21201124509 (SL350)

Range of filters



Magnifications

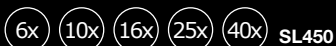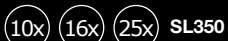

Optometry





WIDE OBSERVATION

- Galilean binocular optical system offering wide field of view and qualitative images.
- Magnification levels easily accessible by drum rotation.



BRIGHT ILLUMINATION

- Tower illumination systems, featuring up to 20° tilting tower illumination.
- Continuous adjustment of slit aperture.
- Available with halogen or LED illumination source.



INTEGRATED BASE

- Easy installation on refraction unit or instrument table thanks to the integrated base and chinrest.

OPTIONS

- Yellow filter, polarised filter, tonometer adaptor.

MAGNIFICATION		SLIT	
Type	Galilean	Slit width	Continuous from 0 to 14mm
Cylinder	12.5x	Slit length	Continuous from 1 to 14mm
Selection	Drum rotation	Aperture	0.2 to 14mm
Field of view	10x (22mm) - 16x (13.5mm) - 25x (8.5mm)	Slit angle	0 to 180°
		Tilt	5, 10, 15, 20°
SL350L/SL350H	6x (34.7mm) - 10x (22mm) 16x (13.5mm) - 25x (8.5mm) 40x (5.5mm)	ILLUMINATION LAMP	
		SL350H/SL450H	6V, 20W halogen lamp
		SL350L/SL450L	3v,3W LED lamp
SL450L/SL450H		FILTERS	
		Heat-absorbing, ND, Red-free, Cobalt blue	
BASE		POWER-SUPPLY	
Longitudinal movement	110 mm	220/110V (±10%)	
Lateral movement	110 mm	DIMENSION AND WEIGHT	
Vertical movement	30 mm	450 (W) x 500 (D) x 740 (H) mm 22kg	
Chinrest (vertical)	80 mm	☑️ Conformity marking	

As improvements are made, these specifications and pictures are not contractually binding and may be changed without prior notice.



Instruments Division
64 bis, avenue Aubert
94306 Vincennes Cedex France
Tel: 33 (0)1 72 94 71 00
Fax: 33 (0)1 72 94 70 27
www.essilor-instruments.com

Distributed by:



PT. Meta Mata Medika www.metamedika.com
MASPION PLAZA, Jl. Gunung Sahari Raya Kav.18, 3rd Floor Unit 3F
Jakarta Utara, INDONESIA 14420
Tel. + 62 21 6470 1108, Fax: +62 21 6470 1114 • Email: info@metamedika.com

Your local distributor: